

PCR Cooler, 96-Well

Ice-Free Cooling System for PCR Preparation and Storage

Cools industry standard format PCR tubes, strips and single format as well as most standard and low profile 96-Well plates

Consistent & Reliable

- Modern design cools faster in freezers and stays cold longer on the bench top
- Maintains the entire 96-Well plate at a temperature of 0 °C for over an hour
- Maintains sample temperature of less than 5°C for up to three hours (with lid off)
- When the temperature ceiling of 6°C has been exceeded, users are signaled by a contrasting color change of the cooler or a single well

Alphanumeric Well Identification
for quick location and retrieval of tubes



Lid doubles as a stand
enables the rack to stay cold longer as there is no conductive heat transfer from surface beneath




Vented Lid
enables faster cooling of samples while in freezer



Ergonomic Textured Finger Grips
provide a secure hold, reducing the chance of dropping the PCR rack, and makes handling and transportation easy, especially when using gloves




Two Color Contrasting Options Available
purple to pink or dark blue to light blue allowing easy monitoring of samples

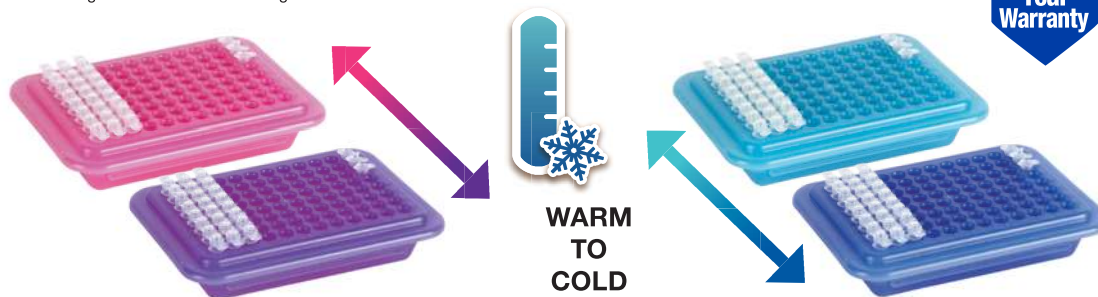

US Patent Pending, Community Design No: 004101756-0001-0002 Pack of 2

Item No.	Warm Color	Cold Color	Wells	Tube Sizes mL	L x W x H in	L x W x H cm
120727	Pink	Purple	96	0.2 Assay Plates	4.3 x 5.8 x 1.9	11 x 14.8 x 4.8
120728	Light Blue	Dark Blue				
120729	Assorted*					

*Assorted colors (Light Blue/Dark Blue & Pink/Purple)

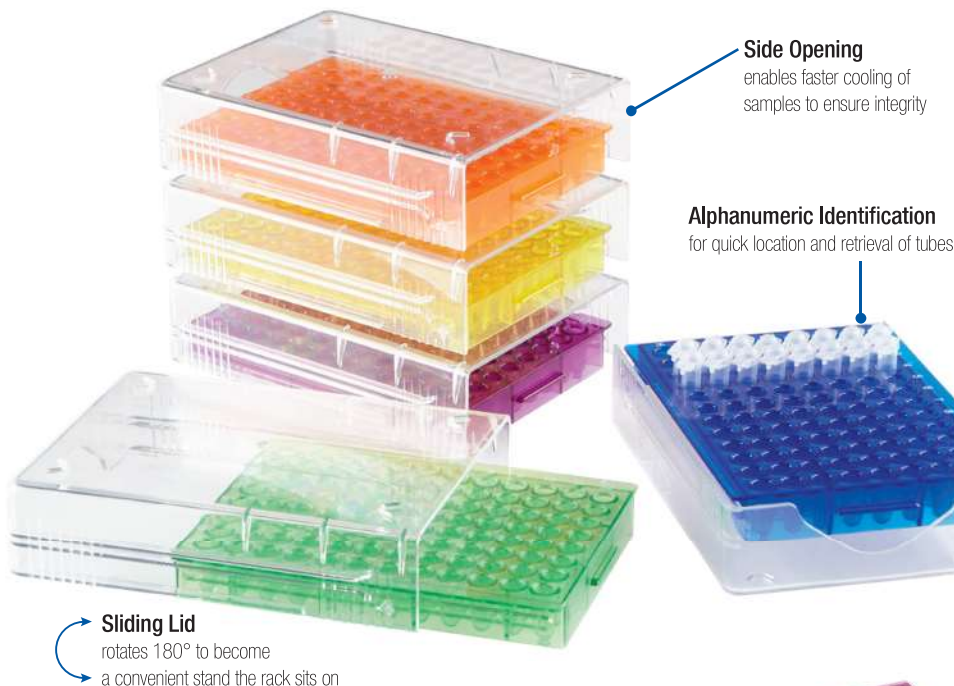

Low-Temp PCR Rack: 96-Well

Slide, set, work & store PCR samples

Multi-purpose rack converts from compact workstation to lidded storage solution accommodating individual and strip 0.2 mL tubes

Versatile & Robust Solution

- Polycarbonate lid serves as cover and elevated rack base - just slide open and pivot lid under rack for compact workstation
- Lid can be used in right or left-handed orientation and is detachable
- Stackable saving valuable bench-top and freezer space
- Multiple colors offered for coding options



US Patent Pending, Community Design No: 004101756-0001-0002 Pack of 5

Item No.	Color	Item No.	Color
120538	Assorted*	120540	Green
120539	Blue	120541	Purple

*Assorted colors (Purple, Blue, Green, Orange, Yellow)

Tubes mL	Wells	Well Dia. mm	L x W x H in	L x W x H cm
0.2	96	5.2	4.9 x 3.6 x 1.5	12.5 x 9.2 x 3.7